

UltraEco ULT Freezer

Redefining Cold Storage Excellence



World's No.1
Energy Efficiency



Ultra-Low
Carbon Footprint



Superior
Performance



Uncompromising
Safeguards



User-Centric
Convenience

*Based on ENERGY STAR certified ULT freezer energy data (as of Feb. 2026)



DW-86L580BPT/730BPT/ 830BPT/960BPT

Qingdao Haier Biomedical Co., Ltd.

No.280 Feng Yuan Road, High-tech Zone,
Qingdao, 266111, P.R. China
E-mail: inquiry@haierbiomedical.com
Website: www.haiermedical.com

Apr, 2026



Haier Biomedical



Haier Biomedical
International



Haier Biomedical
International



Haier Biomedical
International



Haier Biomedical
International


Optimal Energy Efficient Performance



The UltraEco ULT freezer – an ultra-low impact, greener alternative for cold storage

The UltraEco ULT freezer combines a hydrocarbon refrigerant, inverter compressor technology and new sustainable material science to significantly reduce electricity consumption:

Power Consumption (kWh/year/2"cryobox)

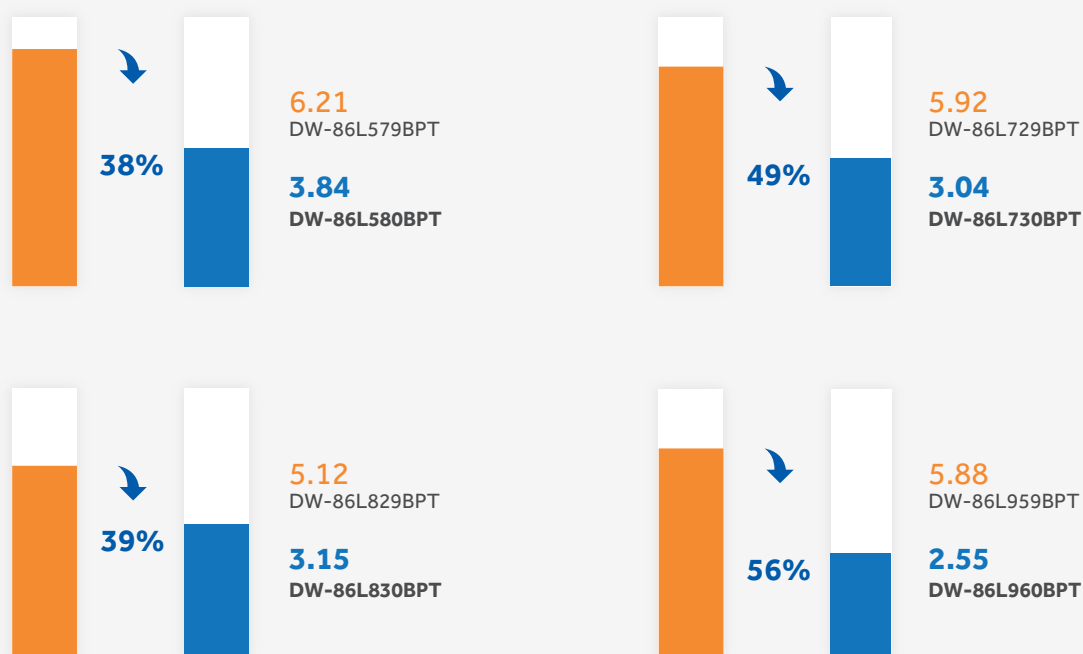

Enhanced Performance



An upgraded compressor, enhanced refrigeration system and new sustainable insulation material deliver a range of benefits compared to the previous generation of ULT freezers.

↓ 49%
Decrease in Energy Consumption

🕒 80%
Faster Re-Entry Time

🔊 42dB(A)
Lower Sound Level

🕒 15 min
1 Min Door Opening Recovery Time

🕒 1 hour
Improvement in warm-up time from -80 to -50°C

🌡 ±2°C
Outstanding Temperature Uniformity

*Above data based on DW-86L730BPT

Green and Environmental Friendly



Naturally occurring hydrocarbon refrigerants and insulation materials are a feature of the UltraEco ULT freezer – the first of its kind in this sector. These non-toxic, eco-friendly materials are more sustainable, non-ozone depleting and reduce the carbon footprint, providing a greener, more environmentally- friendly insulation solution than other current generation freezers.

Safety and Reliability

- Dual sensors deliver precise temperature change detection, doubling sample storage security.
- The compressor's self-diagnostic function triggers instant alarms upon malfunctions, enabling prompt maintenance to prevent sample damage.
- Optimized VIP insulation ensures extended thermal retention, safeguarding samples during brief power outages.
- Optional built-in backup system and independent control panel guarantee sample safety, even during extended power outages.

Ergonomic Design



10.1 Inch Capacitive Touch Screen

State-of-art UI design coupled with sample management system to provide the optimal user experience



Innovative Light Bar Design

Clearly indicate the operating status of the freezer and offer strengthened visual alarm indications for effortless at-a-glance monitoring



Improved Pressure Relief Port

Re-entry time shortened to just 60 seconds







USB Interface

Enables users to download historical temperature data for compliance/auditing purposes


Performance Curve

Temperature Pull-down & Pull-up Curves at 20°C Ambient

— Pull-down
— Pull-up

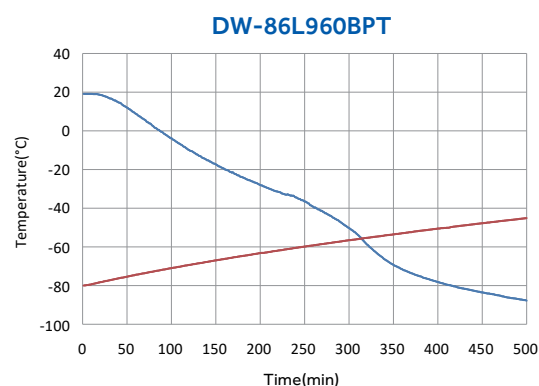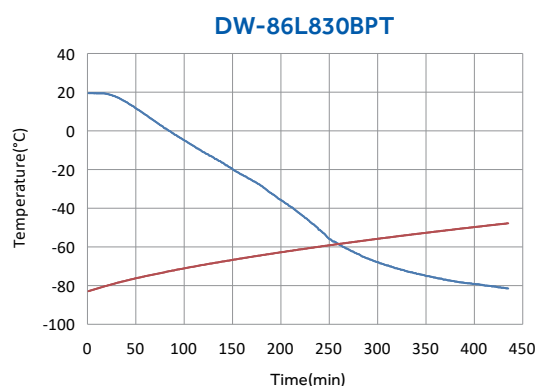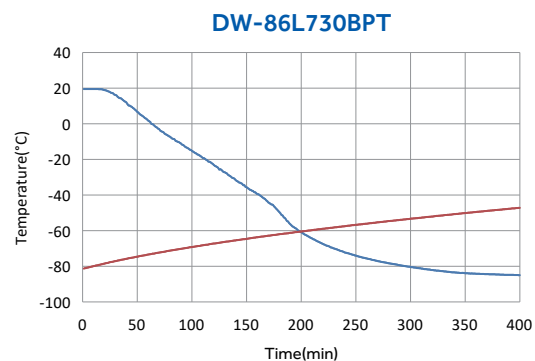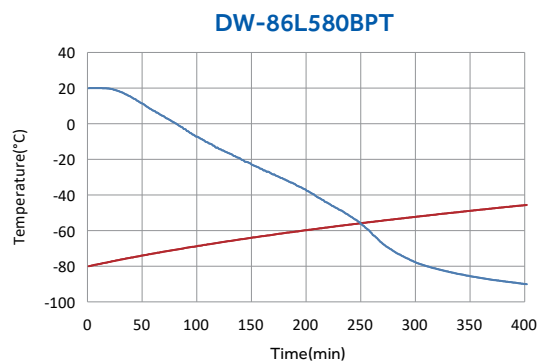

Product Dimensions

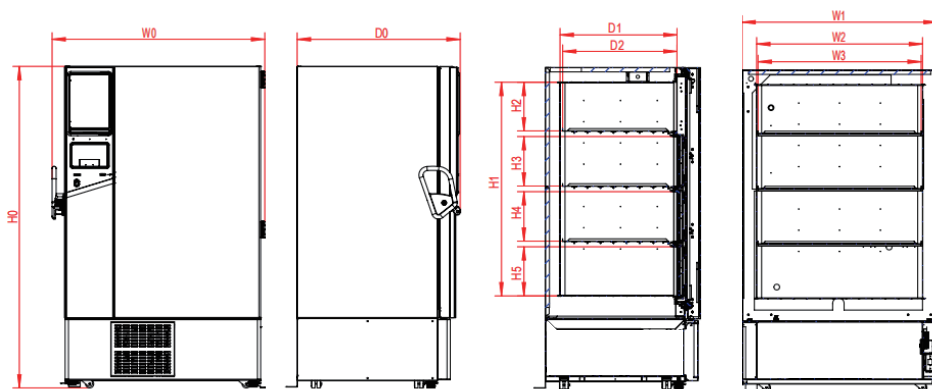

Remark

- W0 ----- Overall Unit Width
- W1 ----- Cabinet Exterior Width
- W2 ----- Inner Liner Width
- W3 ----- Shelf Clear Width (for freezer racks)
- D0 ----- Overall Unit Depth
- D1 ----- Inner Liner Depth
- D2 ----- Maximum Item Storage Depth
- H0 ----- Overall Unit Height
- H1 ----- Inner Liner Height
- H2 ----- 1st Layer Height
- H3 ----- 2nd Layer Height
- H4 ----- 3rd Layer Height
- H5 ----- 4th Layer Height

Code Model	W0/mm	W1/mm	W2/mm	W3/mm	H0/mm	H1/mm	H2/mm	H3/mm	H4/mm	H5/mm	D0/mm	D1/mm	D2/mm
DW-86L580BPT	895	800	620	604	1970	1310	298	303	303	301	998	716	700
DW-86L730BPT	1046	946	766	750	1970	1310	298	303	303	301	998	716	700
DW-86L830BPT	1145	1050	870	854	1970	1310	298	303	303	301	998	716	700
DW-86L960BPT	1296	1196	1016	999	1970	1310	298	303	303	301	998	716	700

Specifications

Model			DW-86L580BPT	DW-86L730BPT	DW-86L830BPT	DW-86L960BPT
Technical Data	Cabinet Type		Upright	Upright	Upright	Upright
	Ambient Temperature (°C)		10-32	10-32	10-32	10-32
	Cooling Type		Direct Cooling	Direct Cooling	Direct Cooling	Direct Cooling
	Defrost Mode		Manual	Manual	Manual	Manual
	Refrigerant		HC	HC	HC	HC
	Sound Level (dB(A))		≤42	≤42	≤46	≤40
Performance	Cooling Performance (°C)		-86	-86	-86	-86
	Temperature Range (°C)		-40~-86	-40~-86	-40~-86	-40~-86
Control	Controller		Microprocessor	Microprocessor	Microprocessor	Microprocessor
	Display		LCD Touchscreen	LCD Touchscreen	LCD Touchscreen	LCD Touchscreen
Electrical Data	Power Supply (V/Hz)		100-230/50/60	100-230/50/60	100-230/50/60	100-230/50/60
	Max Power (W)		800	900	1050	900
Dimensions	Capacity (L/Cu.Ft)		580/20.5	730/25.8	830/29.3	960/33.9
	Net/Gross Weight	kg	287/326	315/368	336/390	375/426
		lbs	632.7/718.7	694.5/811.3	740.8/859.8	826.7/939.2
	Interior Dimension (W×D×H)	mm	620×716×1310	766×716×1310	870×716×1310	1016×716×1310
		in	24.4×28.2×51.6	30.2×28.2×51.6	34.3×28.2×51.6	40.0×28.2×51.6
	Exterior Dimension (W×D×H)	mm	895×998×1970	1046×998×1970	1145×998×1970	1296×998×1970
		in	35.2×39.3×77.6	41.2×39.3×77.6	45.1×39.3×77.6	51.0×39.3×77.6
	Packing Dimension (W×D×H)	mm	960×1150×2160	1110×1150×2160	1220×1150×2160	1400×1150×2160
		in	37.8×45.3×85	43.7×45.3×85	48×45.3×85	55.1×45.3×85
	Cabinet Width (without handle and hinge)		mm 800	946	1050	1196
Functions	Container Load (20'/40'/40'H)		12/24/24	10/20/20	8/20/20	8/16/16
	High/Low Temperature Alarm		Y	Y	Y	Y
	Hot Condenser Alarm		Y	Y	Y	Y
	Power Failure Alarm		Y	Y	Y	Y
	High/Low Voltage Alarm		Y	Y	Y	Y
	Sensor Error Alarm		Y	Y	Y	Y
	Low Battery Alarm		Y	Y	Y	Y
	High Ambient Temperature Alarm		Y	Y	Y	Y
	Door Ajar Alarm		Y	Y	Y	Y
Accessories	Start-up Delay		Y	Y	Y	Y
	Caster		Y(4)	Y(4)	Y(4)	Y(4)
	Leveling Feet		Y(2)	Y(2)	Y(2)	Y(2)
	Porthole		2×rear, (15mm)diameter	2×rear, (15mm)diameter	2×rear, (15mm)diameter	2×rear, (15mm)diameter
	Inner Doors		4	4	4	4
	Shelves		3- stainless steel adjustable shelves	3- stainless steel adjustable shelves	3- stainless steel adjustable shelves	3- stainless steel adjustable shelves
	Shelf Dimensions (W×D)(mm)		606×650	752×650	856×650	1002×650
	Max. Load per Shelf (kg)		72	90	108	126
	Freezer Racks Qty (2" boxes)		16	20	24	28
	Capacity (2" boxes)		400	500	600	700
	USB Interface		Y	Y	Y	Y
	Remote Alarm		Y	Y	Y	Y
	5V Power Supply Port		Y	Y	Y	Y
	Temperature Chart Recorder		Optional	Optional	Optional	Optional
	RS232/485 Port		N/A/Standard	N/A/Standard	N/A/Standard	N/A/Standard
	CO ₂ Backup System		Optional	Optional	Optional	Optional
	LN ₂ Backup System		Optional	Optional	Optional	Optional
Other	Certifications		CE, UL, ENERGY STAR	CE, UL, ENERGY STAR	CE, UL, ENERGY STAR	CE, UL, ENERGY STAR

*Haier Biomedical reserves the right to change products and specifications without prior notice.